

AGELOC® WELLSPA iO™



ageLOC WellSpa iO delivers more of the benefits every body needs to thrive. It's the personalized wellness and beauty device and self-care system to restore, revitalize, and relax. All so you can look and feel like your whole best self.

MEET AGELOC WELLSPA iO

Restore your overall wellness to feel more balanced, revitalize your skin for a self-confidence win, and relax after exercise or a long day with the power of ageLOC WellSpa iO—our personalized wellness and beauty device system. The exclusive design features adaptive microcurrent technology that comfortably adjusts to your skin in real-time, an IoT-connected app to keep you meeting and exceeding your goals, and powerfully formulated ageLOC products that produce immediate results you can feel and see.

Unleash the freedom to feel whole again with ageLOC WellSpa iO.

WHY YOU'LL LOVE IT

RESTORE BENEFITS

Restores overall feeling of wellbeing for better holistic wellness.

- Helps depuff the skin's appearance through skin stimulation.
- Supports overall health and wellbeing.
- Helps you feel more centered and balanced.
- Promotes stimulation and purification to refresh the body's skin.

REVITALIZE BENEFITS

Revitalizes your appearance for better body confidence.

- Visibly contours by minimizing the appearance of cellulite for smoother, radiant skin.
- Visibly firms and tightens skin for a more youthful appearance.
- Helps to nourish the skin to promote smoother and healthier-looking skin.
- Visibly tones and sculpts through skin invigoration.
- Evens skin tone for a visibly brighter, vibrant skin
- Hydrates and helps reduce crepey looking skin.

RELAXING BENEFITS

Relaxes your body for better pre- and post-workout rejuvenation.

- Provides a great pre-workout massage.
- Provides a great post-workout massage.
- Encourages overall skin relaxation and leaves you feeling refreshed.
- Helps your skin feel relaxed.
- Provides a gentle and comfortable massage while contributing to skin relaxation.
- Provides a relaxing massage that leaves a cooling sensation on the skin.
- Provides a smooth, comfortable massage while contributing to skin stimulation.

DEVICE FEATURES

- The surface topography provides a smooth, comfortable experience while contributing to skin stimulation.
- The nodes and current together encourage stimulation.

AGELOC® WELLSPA iO™

- The patent-pending adaptive microcurrent continuously and comfortably adjusts to your skin in real time—for an optimized, just-for-you treatment, every time.
- The silicone sleeve was designed to make cleaning simple and minimize mess.

NU SKIN VERA® APP FEATURES


- Designed to customize and complement as you restore, revitalize, and relax, to help you conquer your everyday wellness and beauty goals.
- Tutorials to guide you every step of the way.
- Delivers complementary product recommendations to help you reach your personal skincare goals.
- Easy reordering so you always have your products on hand.
- Designed to provide a personalized iO Score to help you stay motivated and reach your wellness and beauty goals.
- Developed to track your unique usage for fun, incentivizing badges.
- Ongoing optimization through both app and device updates.

CLINICAL RESULTS THAT PROVE IT REVITALIZE STUDY

Thirty-three female subjects between 25–60 years of age with mild to moderate cellulite on the back of their upper thigh, back of upper arms and/or abdomen were recruited to use the ageLOC WellSpa iO Revitalize Regimen. This regimen consisted of using ageLOC Body Polish as an exfoliating cleanser in the above areas, the WellSpa iO device with ageLOC Body Serum for three minutes, three times per week, and ageLOC Dermatic Effects® twice daily. Each area was evaluated by clinical grading and self-perception at baseline and weeks two, four, eight, and 12. The back of the upper thigh was also evaluated by instrumentation at these same time points.


The Cutometer results illustrating the improvements in skin firmness and elasticity on the thigh are presented below.

Instrumentation: Cutometer

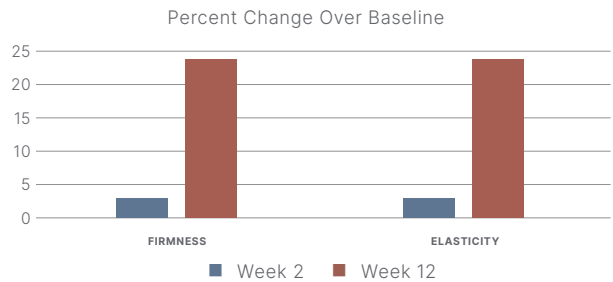

The clinical grader assessed the appearance of cellulite for each area. Below is a representation of the changes seen on the thigh.

Clinical Grading: Appearance of Thigh Cellulite

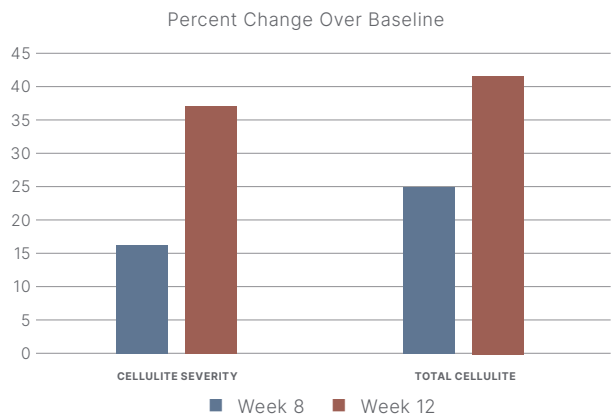

Subjects also noted improvements in their skin's appearance within the following attributes for each area:

- Cellulite
- Firmness
- Contouring
- Tightness
- Smoothness
- Hydration
- Overall appearance

Our research has also illustrated that with increased frequency and timing of using WellSpa iO with ageLOC Body Serum, benefits can be seen more quickly and with more significant results through 12 weeks.

AGELOC® WELLSPA iO™

HOW TO USE IT

Before use, please refer to the ageLOC WellSpa iO manual for safety information and the Vera® app for usage directions.

Whichever routine you choose, we formulated our luxurious, leave-on products with ageLOC technology that targets the sources of aging to preserve the look of youth—even after your treatment is done.

AGELOC WELLSPA IO + AGELOC BODY SERUM

Use three times per week to daily. After cleansing, apply a generous amount of ageLOC Body Serum onto problem area to be treated (e.g., upper arms, upper thigh/buttock, or abdomen) and begin treatment with WellSpa iO, focusing on upward/inward movements, gliding up in slow, circular motions and massaging back to the start. Reapply as needed during the session to ensure smooth, continuous connectivity between the skin and device. Massage any remaining serum into the skin.

For best results, cleanse with ageLOC Body Polish as the first step, followed by ageLOC Dermatic Effects® at the end of your routine.

AGELOC WELLSPA IO + AGELOC BODY ACTIVATING GEL

Apply a generous amount of ageLOC Body Activating Gel to the desired treatment area up to once a day. Focusing on upward/inward movements, glide the device straight upward or in slow, circular massaging movements upward and lightly back to the start. Use with your ageLOC WellSpa iO device. Massage any remaining gel into your skin.

Note: Connect to the Nu Skin Vera app or see the ageLOC WellSpa iO manual for complete usage details.

WHAT POWERS OUR AGELOC BODY PRODUCTS

AGELOC BODY SERUM

ageLOC Body Serum is used with the WellSpa iO device and was formulated to target visible contouring, reduction in cellulite appearance, and improved skin brightness and overall appearance. Not only can it help to visibly firm and tighten the skin, but it also has great hydrating ingredients that help reduce visible skin crepiness, leaving your skin smooth and more youthful looking.

- Body Silhouette Blend—a unique blend of botanical ingredients that targets multiple concerns related to the appearance of skin on the body—firmness, lifting, contouring, and smoothing the appearance of cellulite.
- AI Peptide—a unique natural oligopeptide identified via AI to have key benefits for skin in firming, smoothing, and antioxidant protection.
- Caesalpinia Spinosa Fruit & Kappaphycus Alvarezii Extract—a weightless-feeling film-former that mechanically improves the lifted and firmed appearance of skin.
- Skin Conditioning Blend (sunflower, macadamia, allantoin, and panthenol)—helps to improve the look and condition of skin by reducing flaking and restoring skin suppleness.
- ageLOC Blend—a proprietary ingredient blend that targets the sources of aging to preserve the look of youth.

AGELOC BODY ACTIVATING GEL

ageLOC Body Activating Gel is a refreshing formula used with your Restore and Relax WellSpa iO routines. It's formulated to work with the ageLOC WellSpa iO device to help you look and feel your best.

- ageLOC Blend—is proprietary ingredient blend that targets the sources of aging to preserve the look of youth.
- WellBe Blend—a unique blend of Crithmum Maritimum and Ashwagandha extracts that help target the overall feeling of wellbeing.

YOU MIGHT ALSO LIKE

- ageLOC Body Polish
- ageLOC Dermatic Effects®
- Beauty Focus® Collagen+
- ageLOC Meta™
- Epoch® IceDancer® Invigorating Leg Gel

LEARN MORE ABOUT IT

How do you integrate wellness and beauty in one device?

We treat our bodies like an ecosystem and want holistic solutions that complement our personal health and evolving needs. ageLOC WellSpa iO does just that.

AGELOC® WELLSPA iO™

From the powerful ageLOC formulas to the unique topography, adaptive microcurrent, and connected iO experience, WellSpa iO provides restoring, revitalizing, and relaxing wellness and beauty benefits that help you look and feel your best

ageLOC Body Serum and ageLOC Body Activating Gel have similar key ingredients—how do they provide different benefits?

Ingredients can have different types of benefits that fit different desires. This is, in part, due to differences in extraction methods for extracts. Or, at other times the same ingredient can have different targets. So, while an ingredient may be in one formula for one reason, it may have different benefits in another formula.

What are the three Rs, and how do they impact my wellbeing?

The ageLOC WellSpa iO system is designed to deliver the power of overall wellness, body confidence, and rejuvenation through our two Rs—Restore, and Revitalize.

- **Restore**—Designed to depuff and energize your skin, and feel your best any time of day as ageLOC WellSpa iO restores a feeling of “you-got-this” total-body health.
- **Revitalize**—ageLOC WellSpa iO helps you look your best as it stimulates your skin to visibly firm, tighten, brighten, tone, and hydrate—while helping reduce the appearance of cellulite.
- **Relax**—the system is also designed to provide a gentle and comfortable massage to relax the skin pre and post workout.

What benefits will I enjoy because of Adaptive Microcurrent Technology?

The short answer is that it's designed to help you maximize your wellness and beauty goals. But to be technical, the microcurrent continuously and comfortably adjusts to your skin in real time—for an optimized, just-for-you treatment every time. Plus, when you pair ageLOC WellSpa iO with our luxurious, leave-on products, the ageLOC technology targets the sources of aging to preserve the look of youth—even after your treatment is done.

How is ageOC WellSpa iO System unique?

WellSpa iO is a complete reimagination of the at-home spa experience. It utilizes a more powerful Adaptive Microcurrent that continuously adapts to the changes in your skin conductivity throughout each session to provide you with an optimal and comfortable experience.

WellSpa iO also incorporates one-of-a-kind, patent-pending metallic nodes that help provide a smooth, comfortable massage while contributing to skin stimulation—designed to perfectly fit the contours of your body and glide smoothly over your skin. The accompanying topicals (ageLOC Body Serum and ageLOC Body Activating Gel) are a seamless leave-on experience packed with powerful ingredients. At the same time, Body Polish provides complementary benefits of cleansing and polishing the skin for a smooth, radiant finish. This system not only helps to contour and beautify the skin, but it also revitalizes for an overall uplifting experience as you personalize your wellness and beauty journey.

How does the microcurrent in WellSpa iO differ from ageLOC Galvanic Body Spa®?

ageLOC Galvanic Body Spa utilizes a pulsing direct current to stimulate the skin. WellSpa iO utilizes Adaptive Microcurrent. This is a biphasic pulsing current that continuously adapts to the conductivity of the skin in real time. The ability to adjust to the skin's conductivity allows the device to provide the optimal current intensity without compromising comfort.

What is the maximum amount of time recommended to use WellSpa iO?

You may use your ageLOC WellSpa iO device system up to 5 minutes, per body area, per day.

Why are there different usage movements for WellSpa iO?

Those looking to restore their wellbeing target an entire area (e.g., the whole thigh, the entire arm) and focus on basic upward and inward movements over the limb. When looking to revitalize your appearance, users will want to use movements that are directed up and/or inward, but also provide some gentle massaging of the area to help improve the appearance of the skin. These areas are likely more targeted as opposed to the entire limb.

AGELOC® WELLSPA iO™

Why do the metallic nodes have such a unique shape?

The WellSpa iO device's one-of-a-kind metallic node design was created to help provide a smooth, comfortable massage while contributing to skin stimulation. The device was designed to perfectly fit to the contours of your body and glide smoothly over your skin.

Do I have to connect with the app to use my ageLOC WellSpa iO?

You'll still enjoy all the great wellness and beauty benefits if you don't connect. However, ageLOC WellSpa iO is created to connect so you'll want to benefit from all the features developed exclusively for you, including personalized goal setting, guided tutorials, custom routines, and usage tracking.

How do I know which routine is right for me?

You know your wellness and beauty needs best, and there are plenty of routines to choose from once you're in the Vera® app. Feel free to try them all to get all the amazing benefits we're sure your skin and body will love.

How does ageLOC WellSpa iO support Nu Skin's commitment to sustainability?

As responsible stewards of our planet, we have proudly designed the WellSpa iO with a rechargeable battery and to be more water resistant, increasing the device's longevity. In addition, the carton is made from 100% recycled material. Please check your local device recycling capabilities and recycle where possible. By purchasing WellSpa iO, you're joining us in taking an active stride to ensure the future is bright for both the planet and each of us. Please visit www.nuskin.com/sustainability for more information and details on how to recycle.

WHAT'S IN IT

AGELOC BODY SERUM

Water (Aqua), Glycerin, Propanediol, Pentylene Glycol, PEG-16 Macadamia Glycerides, PEG-10 Sunflower Glycerides, PEG-40 Hydrogenated Castor Oil, Camellia Sinensis Leaf Extract, Oligopeptide-1, Sodium Hyaluronate, Sodium Acetylated Hyaluronate, Kappaphycus Alvarezii Extract, Lotus Corniculatus Flower Extract, Crithmum Maritimum Extract, Caesalpinia Spinosa Fruit Extract, Panthenol,

Arginine, Allantoin, Carthamus Tinctorius (Safflower) Seed Oil, Tocopheryl Acetate, Carbomer, Xanthan Gum, Hydroxypropyl Methylcellulose Stearoxy Ether, Caprylyl Glycol, Hydroxyacetophenone, Ethylhexylglycerin, Sodium Benzoate, Disodium EDTA, Fragrance (Parfum), Sodium Hydroxide, Citric Acid.

AGELOC BODY ACTIVATING GEL

Water (Aqua), Butylene Glycol, Glycerin, Polysorbate 20, Withania Somnifera Root Extract, Crithmum Maritimum Extract, Sodium Hyaluronate, Schizandra Chinensis Fruit Extract, Oryza Sativa (Rice) Extract, Lotus Corniculatus Flower Extract, Allantoin, Bisabolol, Panthenol, Sodium PCA, Aloe Barbadensis Leaf Juice, Carbomer, Fructose, Caprylyl Glycol, Chlorphenesin, Disodium EDTA, Phenoxyethanol, Fragrance (Parfum), Aminomethyl Propanol.



Scan the QR code to connect with Vera®

App Store is a registered trademark of Apple, Inc. and Google Play is a registered trademark of Google LLC